



IndelCheck™ CRISPR/TALEN insertion or deletion detection system

Cat. No. ICPE-050

Cat. No. ICPE-200

Cat. No. TENI-050

Cat. No. TENI-200

Cat. No. TPCR-050

Cat. No. TPCR-200

User Manual

GeneCopoeia, Inc.
9620 Medical Center Drive, #101
Rockville, MD 20850
USA

301-762-0888
866-360-9531

inquiry@genecopoeia.com

www.genecopoeia.com

USER MANUAL

IndelCheck™ CRISPR/TALEN insertion or deletion detection system

I. Introduction	3
II. Contents and storage	4
III. Protocol overview.....	6
IV. Experimental Procedures	7
1. Primer design	7
2. Sample preparation.....	7
3. Target site PCR and product processing	8
4. Denaturation and re-annealing	10
5. Digest with T7 Endonuclease I.....	10
6. Gel analysis	10
7. Sequence verification.....	11
V. Notes and Troubleshooting	12
VI. Appendix.....	14
VII. Limited Use License and Warranty	20

I. Introduction

CRISPR/TALEN introduced double-stranded breaks (DSBs) at desired target sites can be repaired by nonhomologous end joining (NHEJ), which is error prone and typically results in small insertions or deletions (indels) near the DSB. The target region is PCR-amplified, and the PCR products are denatured and re-annealed to allow mismatched DNA to form, such as wildtype/indel mutation mismatches or indel mutation 1/indel mutation 2 mismatches. The re-annealed PCR fragments are incubated with T7 endonuclease I, which recognizes and cleaves mismatched DNA. If two shorter bands of the predicted size are generated, it usually means that CRISPR/TALEN has successfully introduced indel mutations at the targeted chromosomal site.

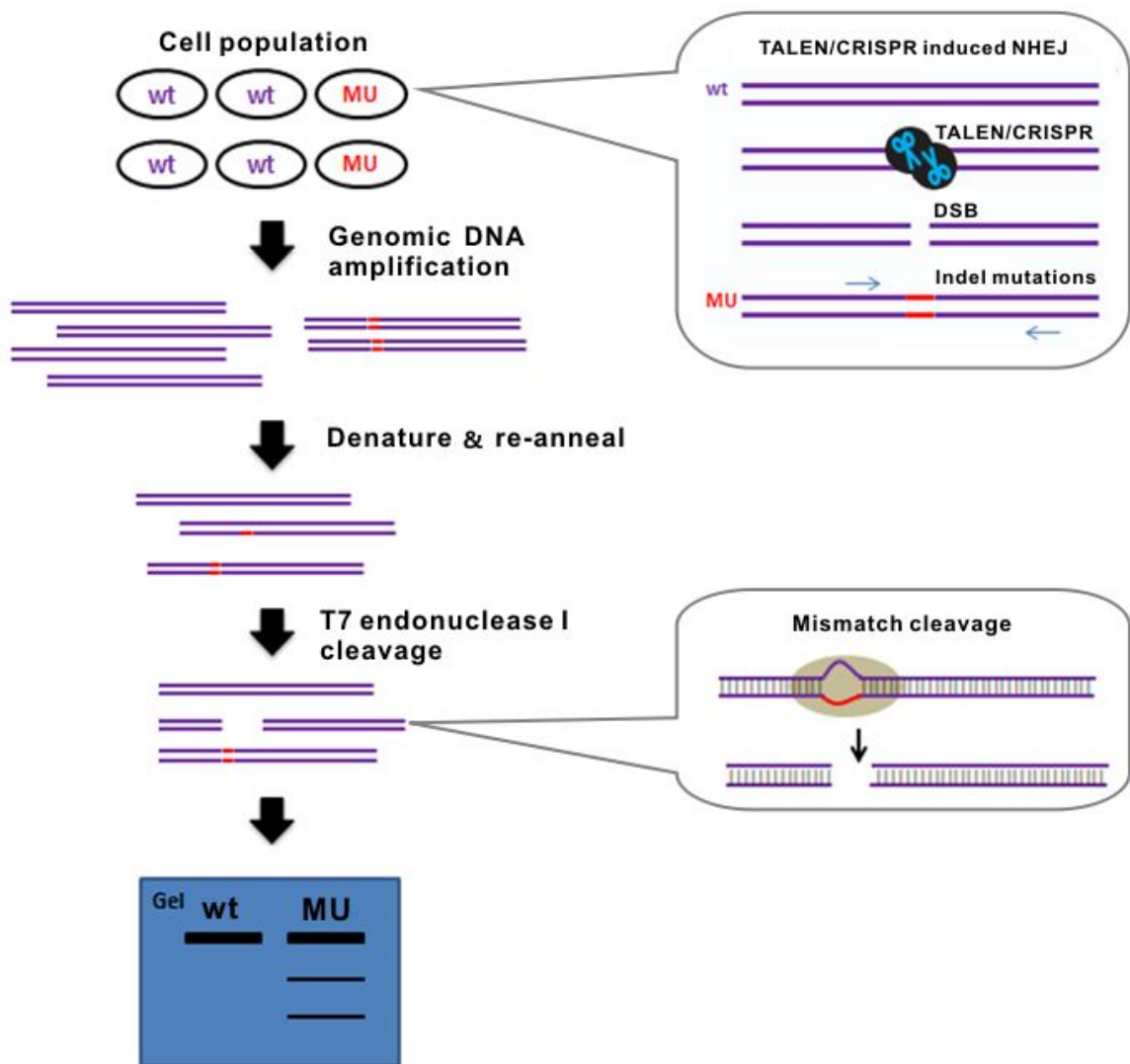


Figure 1. CRISPR or TALEN functional validation using the mismatch cleavage assay

IndelCheck™ CRISPR/TALEN insertion or deletion detection system

The IndelCheck™ CRISPR/TALEN insertion or deletion detection system (ICPE-050, ICPE-200) is a complete system designed to simplify your validation or screening workflow. It includes both the Target Site PCR kit (TPCR-050, TPCR-200) and the T7 Endonuclease I Assay kit (TENI-050, TENI-200). All components are optimized together for the best performance.

The Target Site PCR kit is optimized for robust amplification of the target site. With the proprietary lysis buffer included, genomic DNA isolation is no longer required. The kit contains the high efficiency and high fidelity SuperHeRo DNA polymerase, which produces both blunt-end and sticky-end PCR products compatible with a variety of sequencing vectors. GeneCopoeia also provides the target site PCR primer design and synthesis. Go to our [website](#) for more details.

The T7 Endonuclease I Assay kit contains T7 endonuclease I, which detects and cleaves heteroduplex DNA. It also contains positive controls for both target PCR and indel mismatch cleavage.

Advantages

- Complete system to simplify your CRISPR/TALEN validation and knockout clone screening
- Robust amplification for the target site PCR. No genomic DNA isolation is required
- Easy to use T7 endonuclease I assay with optimized conditions and positive control

II. Contents and Storage

IndelCheck™ CRISPR/TALEN insertion or deletion detection system (ICPE-050, ICPE-200)

The T7 endonuclease I assay kit comes in two sizes:

50-Reaction Kit (Catalog No. TENI-050)

200-Reaction Kit (Catalog No. TENI-200)

Content	Quantity 50 reactions 200 reactions	Shipping temperature	Storage temperature
Digestion reagents			
T7 Endonuclease I (2U/μL)	100μL 400μL	Ice pack	-20°C Stable for at least 12 months
10× T7EN Buffer	100μL 400μL	Ice pack	-20°C Stable for at least 12 months
Control reagents			
Control template & primer mix	200μL 800μL	Ice pack	-20°C Stable for at least 12 months

The target site PCR kit comes in two sizes:

50-Reaction Kit (Catalog No. TPCR-050)

200-Reaction Kit (Catalog No. TPCR-200)

Content	Quantity 50 reactions 200 reactions	Shipping temperature	Storage temperature
Lysis Buffer NEW	650 μ L \times 2 650 μ L \times 8	Ice pack	-20°C Stable for at least 12 months
5XPCR Buffer	250 μ L 1 mL	Ice pack	-20°C Stable for at least 12 months
5XEnhancer	250 μ L 1 mL	Ice pack	-20°C Stable for at least 12 months
20mM Mg ²⁺	125 μ L 500 μ L	Ice pack	-20°C Stable for at least 12 months
25mM dNTP	10 μ L 40 μ L	Ice pack	-20°C Stable for at least 12 months
SuperHeRo DNA polymerase	12.5 μ L 50 μ L	Ice pack	-20°C Stable for at least 12 months

Important note: .

Store all components at -20°C upon receipt.

Materials required but not supplied

The following materials are required but not supplied:


ddH₂O. Avoid using autoclaved H₂O. The recycled steam in some autoclaves can introduce contaminants that may interfere with PCR.

Target site PCR primers specific to target site(s) with T_m \geq 60°C The primers should flank the TALEN or sgRNA target site and generate an amplicon of approximately 500-800 bp with the TALEN or sgRNA target site offset from the center by approximately 100 bp. Make sure the primers are specific for the intended site. Avoid using primers that contain inosine. Avoid to have potential SNP sites or sequence differences between alleles in your target region.

***NOTE:** We highly recommend you to obtain sequence information of the target site of you cell line before primer design and any other experiment. It is possible for mammalian cells, which are usually diploid, to have sequence differences between alleles. This may cause false positive results when using T7 endonuclease I to digest negative controls. Such sites should be avoided when designing target PCR primers. See Appendix 4. for further instructions on checking potential SNP sites using online tools.*

*GeneCopoeia also provides design and synthesis services for sequence-specific target site PCR primers.

IV. Protocol Overview

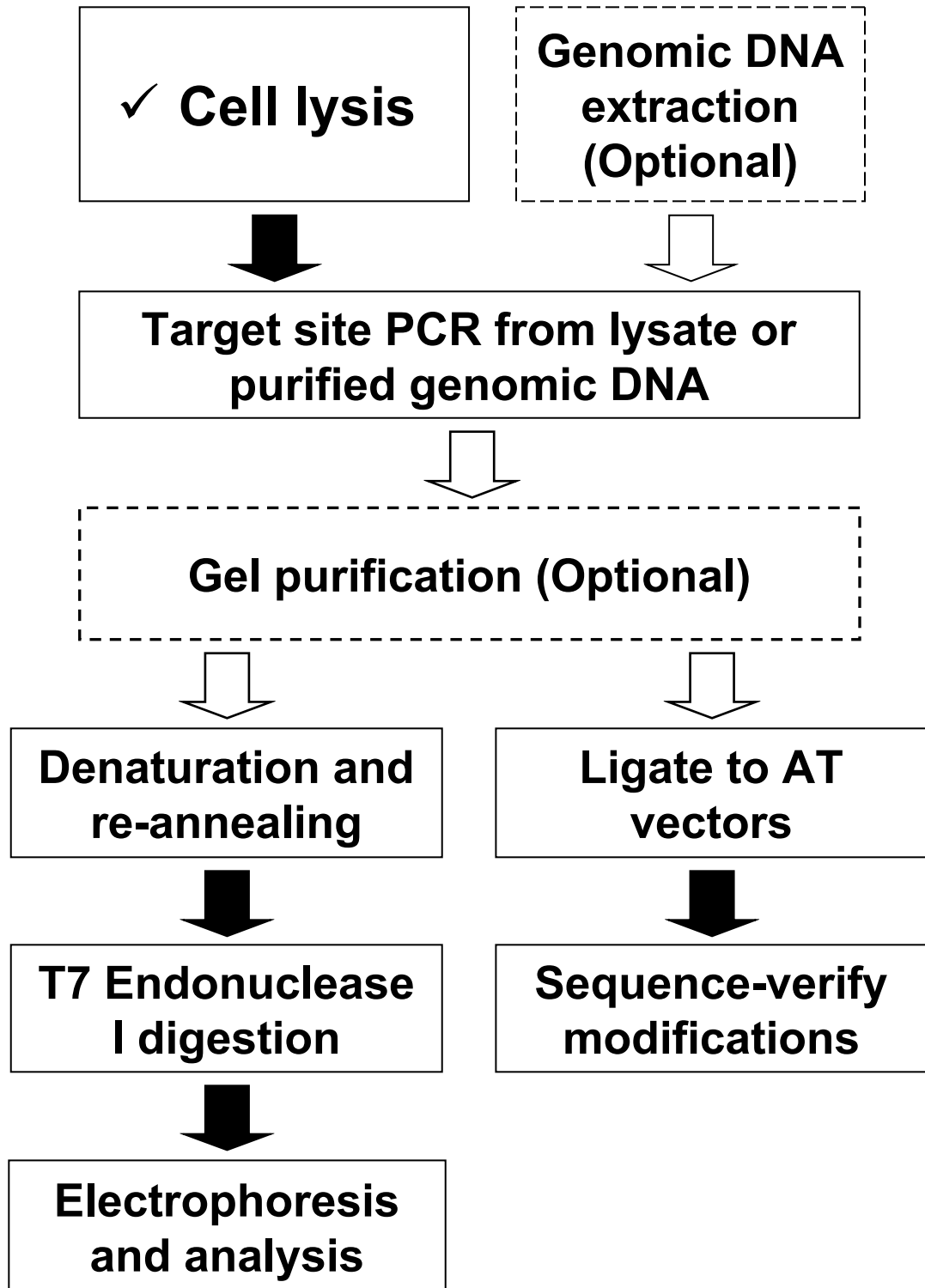

IV. Experimental Procedures

This section provides instructions for validating CRISPR sgRNA or TALEN chromosomal cleavage activity using the IndelCheck™ CRISPR/TALEN insertion or deletion detection system.

In addition to the steps shown in the protocol overview, we also provide brief instructions for PCR primer design and genomic DNA extraction. If you stop before completing all the steps, please store your PCR or digestion products at -20°C until the next step. Avoid repeated freeze-thaw cycles.

We recommend that first time users perform a positive control PCR reaction using the control reagents in the kit. This PCR product can serve as a control for denaturation and re-annealing, as well as mismatch digestion steps.

1. Primer Design

- 1) Target site PCR primers should have a T_m value of no less than 62°C
- 2) For optimum results, the amplicon size range should be approximately 500~800bp.
- 3) Design primers so that the CRISPR sgRNA or TALEN target site is offset from the center of the PCR product by approximately 100 bp. This ensures being able to readily resolve the cleavage products on the gel.

2. Sample Preparation

- Option 1: Genomic DNA extraction
 - a) Harvest cells (no fewer than ~10⁶ cells per well).
 - b) Extract genomic DNA using your method of choice or following the provided protocol of extraction kit manufacturer. Make sure the concentration of genomic DNA solution is above 25 ng/μL.
- Option 2: Cell lysate preparation
 - a) Collect cells from cell culture dish, 6- well plate or 96-well plate. Centrifuge at 3000 rpm at 4°C for 5 min and carefully remove the supernatant.
 - b) Add 300μL 1×PBS. Pipette gently to suspend cells. Centrifuge at 3000 rpm at 4°C for 5 min and remove the supernatant.
 - c) Add 300μL 1×PBS and resuspend cells. Sample the suspension to calculate the cell number if necessary. Centrifuge at 3000 rpm at 4°C for 5 min and remove the supernatant as completely as possible. Proceed to lyse or store the pellet at -80°C.
 - d) Add 25μL Lysis Buffer and lyse cells at 65°C for 15 min, then 95°C for 10 min. Quickly put it on ice afterward.

NOTE: The volume of Lysis Buffer can be adjusted basing on the cell number. At least 50,000 and no more than 5×10^5 cells are recommended for use in 25 μ L Lysis Buffer. For confluent cells of a well of 6-well plate, add 200 μ L~600 μ L Lysis Buffer. For confluent cells of a well of 96-well plate, add 50-100 μ L Lysis Buffer. For amplifying fragment > 1 kb, we suggest prolonging the cell lysis at 65°C for 40 min, but no more than 1h.

However, It is not necessary to obtain complete cell lysis in most experiments. The remaining cells can be stored at -80 °C or for continued culture.

- e) Frozen centrifuge at 12000 rpm for 1 min.

NOTE: Too much floc after centrifuge suggests too little lysis occurred. After transferring the suspension to another tube, the precipitate can be resuspended by adding another 25 μ L Lysis Buffer.

- f) Proceeded to PCR reaction with Target PCR kit. The cell lysate can be stored at 4 °C for no more than one week or -20 °C for several months until use.

3. Target PCR and product processing

1) Target PCR

- a) For PCR from extracted genomic DNA, prepare a Master Mix with target PCR primers flanking the insert as follows:

Item	Amount
5X PCR Buffer	5 μ L
5X Enhancer (optional)	5 μ L
genomic DNA	50-200 ng
Primers(2 μ M)	2 μ L
25mM dNTP	0.2 μ L
SuperHeRo DNAPolymerase (5U/ μ L)	0.25 μ L
20mM Mg ²⁺	2.5 μ L
ddH ₂ O	to 25 μ L
Final	25 μL

- b) For PCR from lysate, prepare a Master Mix with target PCR primers as follows:

Item	Amount
5X PCR Buffer	5 μ L
5X Enhancer (optional)	5 μ L
lysate	1* μ L
Forward and reverse primers (2 μ M)	2 μ L
25mM dNTP	0.2 μ L
SuperHeRo DNApolymerase (5U/ μ L)	0.25 μ L
20mM Mg ²⁺	2.5 μ L
ddH ₂ O	to 25 μ L
Final	25 μL

**To avoid insufficient PCR amplification, adjust the volume of cell lysate based on the cell number. to ensure that no less than 2000 copies of template are present in the reaction. For example, for HT1080 cells, which contains 2 copies of each chromosome, at least 1000 lysed cells are needed in a PCR reaction. To get a bright and clear band on agarose gel, about 10,000 lysed cells are needed.*

NOTE: See Appendix 3. for PCR system for the control mix.

- c) Proceed with PCR using the following program

94°C	5 min	1 cycle
94°C	30 s	35 cycs
58°C	30 s	
72°C	1 min	
72°C	5 min	1 cycle

NOTE: PCR should produce a sufficiently high yield of a SINGLE amplified band of the correct size. We strongly recommend use of a high fidelity DNA polymerase to reduce the amount of base misincorporation during PCR and subsequent false positives. The PCR product can be directly used as the substrate for T7 Endonuclease I digestion.

- 2) Purification or gel extraction of correct-sized band from non-specific PCR background. See Appendix 2 for using gel purification to optimize cleavage of non-specific amplicons

4. Denaturation and re-annealing

- 1) For purified genomic PCR product:

DNA (>25ng/μL)	200~500 ng
10X T7EN Buffer	2 μL
Nuclease-free water	Add up to 19 μL
Total	19 μL

For PCR products amplified with GeneCopia target PCR kit, combine the following

Unpurified PCR product	200~500 ng(5~19μL)
Nuclease-free water	Add up to 19 μL
Total	19 μL


- 2) Mix and centrifuge for a few seconds.
- 3) Heat at 95°C for 5 minutes.
- 4) Reanneal by allowing the denatured PCR products to cool down to RT.

5. Digest with T7 Endonuclease I

- 1) Add 1μL of 2U/μL T7 Endonuclease I.
- 2) Incubate at 37°C for 20-60 minutes.

6. Gel analysis

- 1) Add 1/10 volume of 10 × loading buffer [with 0.1% SDS] to each reaction and mix thoroughly.
- 2) Load half of the mixture into the well of 2% agarose/EtBr gel and run in TAE or TBE buffer.
- 3) Also load a 100bp DNA ladder (#M01010A) as a size reference marker in one of the adjacent wells. Run the gel at 5 V/cm – 11 V/cm until the bromophenol blue reach 2/3 of the length of the gel.

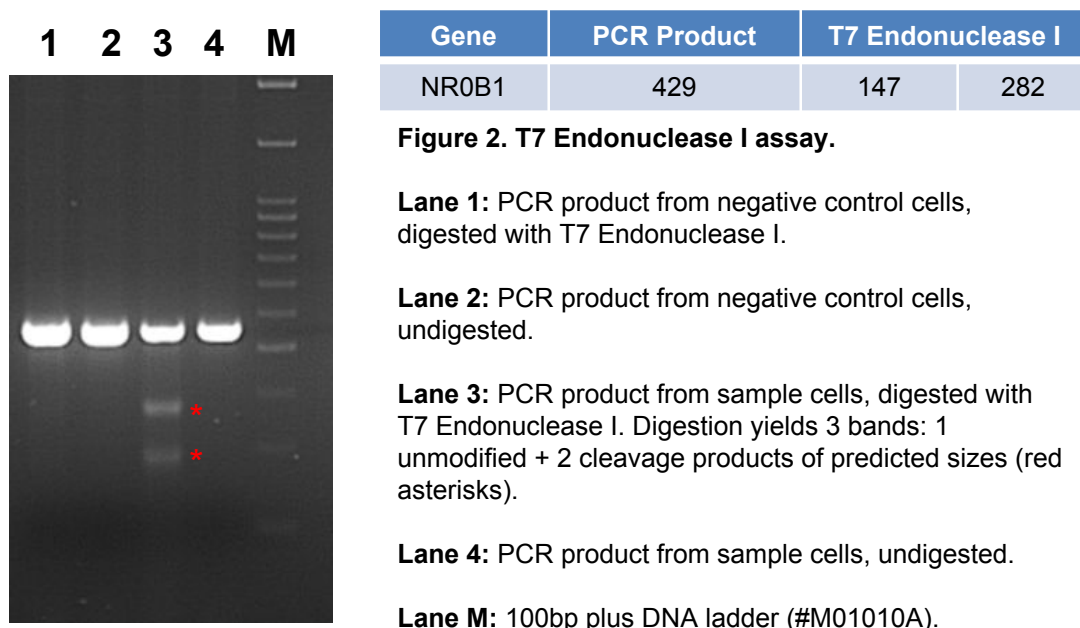

***NOTE:** See Appendix 3. for T7 endonuclease I assay of control.*

7. Sequence verification (following step 3)

- 1) Ligate the digested AT vector backbone (or digested blunt-end vector backbone) and PCR product according to the guidelines provided by the manufacturer of the ligation system. Include a vector-only ligation control..
- 2) Transform the ligated plasmid into a competent E. coli strain, according to the protocol supplied with the cells. Incubate the mixture overnight at 37 °C.
- 3) Pick a minimum of 48 clones to inoculate in 3-5 ml of LB-ampicillin.
- 4) Isolate the plasmid DNA from overnight cultures by using a commercial miniprep kit.
- 5) Verify the sequence of each colony by sequencing from the AT vector backbone using the relevant primers. Reference the sequencing results against the expected genomic sequence to check for the presence of CRISPR sgRNA or TALEN-induced NHEJ modifications. Calculate the percentage of editing efficiency as (no. of modified clones)/(no. of total clones). It is important to pick a reasonable number of clones (>24) to generate an accurate approximation of modification efficiencies.

VI. Notes and Troubleshooting

Cleavage troubleshooting


Problem	Possible Causes	Recommended Solution
Non-specific cleavage bands are observed	Non-specific PCR amplification	<ul style="list-style-type: none"> Do gel purification to ensure that your amplicons are single bands (See Figure 4 in Appendix). Optimize PCR primers on untreated genomic DNA or cell lysates. Optimize PCR condition.
No cleavage products observed	Low T7 Endonuclease I activity	<ul style="list-style-type: none"> If no cleavage is observed for all samples including the positive control, add $MnCl_2$ at a final concentration of 10mM to enhance T7 Endonuclease I activity.
Nonspecific nuclease activity is observed	Incorrect reaction temperature	<ul style="list-style-type: none"> Ensure that samples are incubated at 37°C.
	Reaction time too long	<ul style="list-style-type: none"> Avoid treating with T7 Endonuclease I more than 1.5 hour.
	Poor annealing operation	<ul style="list-style-type: none"> Perform denaturalization and annealing in heated water. Let the reaction cool down naturally with the water. Perform the denaturalization and annealing step in a PCR machine. as follows: (1)95°C 5min (2)94°C(-2°C/cycle), 10-20 sec (3)93°C(-2°C/cycle), 10-20sec and go to step (2),34 cycles
	PCR is introducing mutations	<ul style="list-style-type: none"> Be sure to use a high-fidelity polymerase for PCR amplification
DNA bands are too weak to be observed	Low positive rate of modification	<ul style="list-style-type: none"> If possible, optimize the conditions of your genome editing experiment (e.g., design new CRISPR sgRNA or TALENs).
	Insufficient amount of DNA loaded on gel	<ul style="list-style-type: none"> Be sure to load enough DNA to enable ready visualization on the gel. Also, load equal amounts of total PCR product DNA in each lane.

Target site PCR troubleshooting


Problem	Possible Causes	Recommended Solution
No expected bands	Suboptimal PCR conditions	<ul style="list-style-type: none"> Analyze the sequence of your target. If the GC content is lower than 40%, do not add Enhancer in the PCR reaction.
	Concentration of PCR template is too low	<ul style="list-style-type: none"> Extract and purify genomic DNA to better control the template concentration.
	Incomplete lysis	<ul style="list-style-type: none"> Remove PBS as completely as possible before adding Lysis Buffer, or the remaining PBS will dilute the lysis buffer. Calculate the cell number by hemocytometer or cellometer before lysing the cells. Adjust the volume of lysis buffer according to the cell number.
	Poor PCR primer design	<ul style="list-style-type: none"> Check if the primers match with the template correctly.
Non-specific bands	Suboptimal PCR conditions	<ul style="list-style-type: none"> Increase the annealing temperature to 0~5°C above the T_m value.
	Poor PCR primer design	<ul style="list-style-type: none"> Check the design of PCR primers for possible non-specific binding. If necessary, redesign primers to improve specificity.
	Too much polymerase	<ul style="list-style-type: none"> Decrease the volume of polymerase to 0.2 uL (1U)
Smear	Concentration of PCR template is too high	<ul style="list-style-type: none"> Dilute template 2 or more fold and repeat the PCR.
	Too much polymerase	<ul style="list-style-type: none"> Decrease the volume of polymerase to 0.2 uL (1U)

VII. Appendix

1. Example of using IndelCheck™ CRISPR/TALEN indel detection system to validate CRISPR sgRNA or TALEN cleavage activity

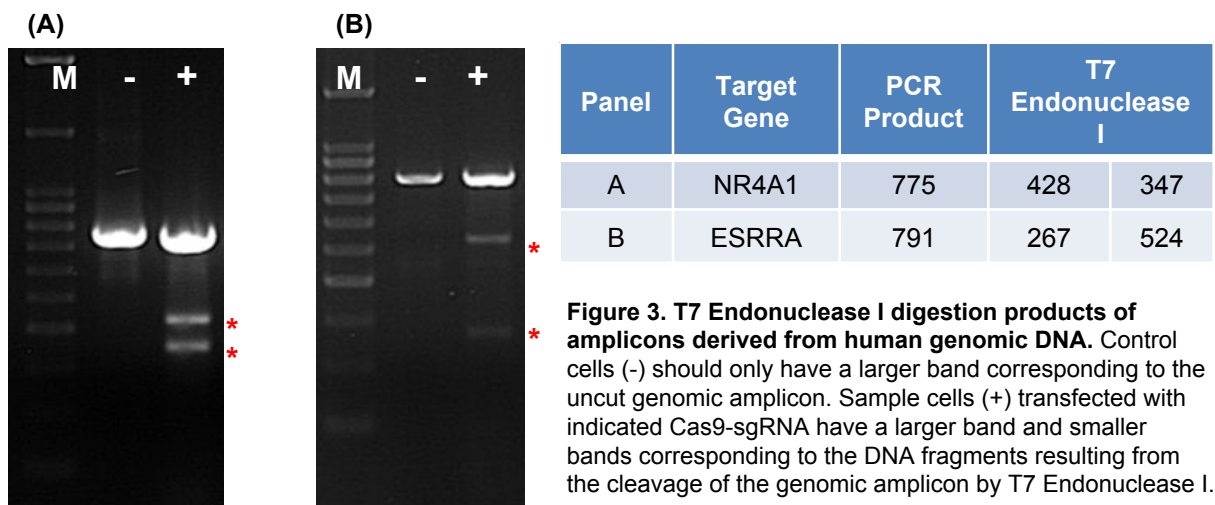

2. Example of using gel purification to optimize cleavage of non-specific amplicons

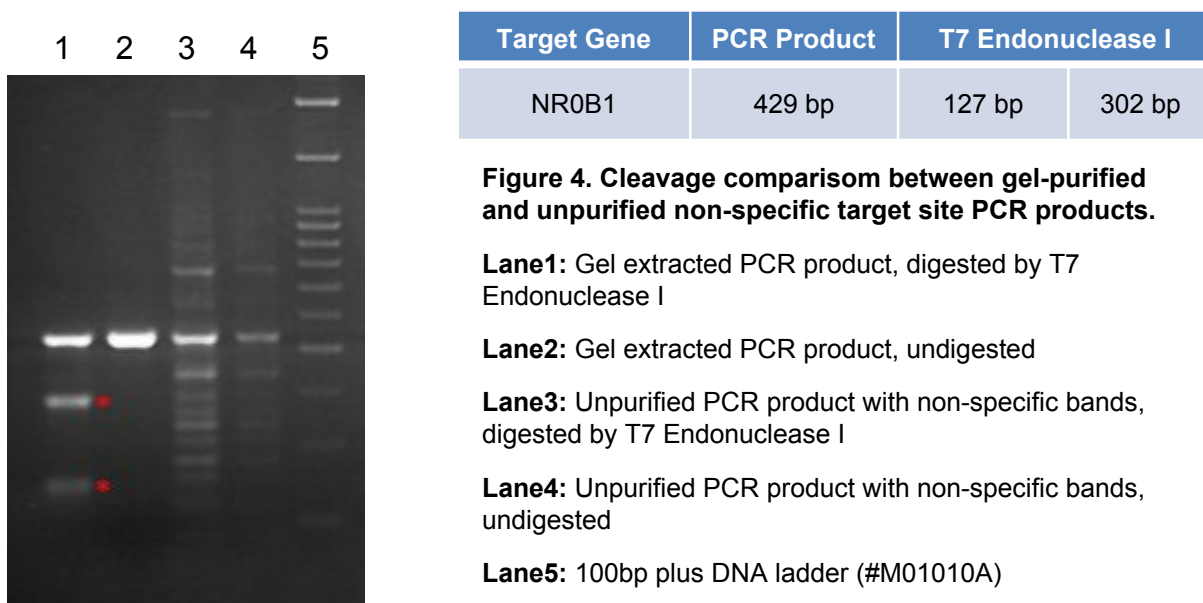

3. Procedures for control target PCR and T7 endonuclease I assay using the control template & primer mix

1) Control PCR

a) Prepare control PCR using the following system.

Item	Amount
5 × PCR Buffer	5 μL
5 × Enhancer (optional)	5 μL
Control template & primer mix	4 μL
25mM dNTP	0.2 μL
SuperHeRo DNAPolymerase (5U/ul)	0.25 μL
20mM Mg2 ⁺	2.5 μL
ddH2O	to 25 μL
Final	25 μL

b) Proceed with PCR using the following program.

94°C	5 min	1 cycle
94°C	30 s	35 cycs
58°C	30 s	
72°C	1 min	
72°C	5 min	1 cycle

2) Denaturation and re-annealing

a) For PCR products amplified with GeneCopeia target PCR kit, combine the following

Unpurified PCR product	200~500 ng(5~19μL)
Nuclease-free water	Add up to 19 μL
Total	19 μL


- b) Mix and centrifuge for a few seconds.
- c) Heat at 95°C 5 minutes.
- d) Reanneal by allowing the denatured PCR products cool down to RT.

3) Digest with T7 Endonuclease I

- a) Add 1µL of 2U/µL T7 Endonuclease I.
- b) Incubate at 37°C for 20-60 minutes.

4) Gel analysis

- a) Add 1/10 volume of 10 × loading buffer [with 0.1% SDS] to each reaction and mix thoroughly.
- b) Load half of the mixture into the well of 2% agarose/EtBr gel and run in TAE or TBE buffer.
- c) Also load a 100bp DNA ladder (#M01010A) as a size reference marker in one of the adjacent wells. Run the gel at 5 V/cm until the bromophenol blue has migrated 2/3 of the length of the gel.

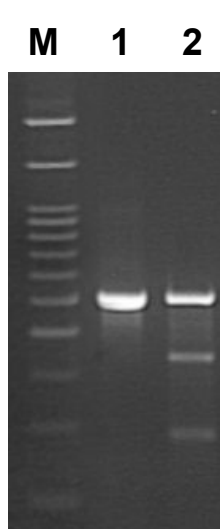

Gene	PCR Product	T7EI	
control	520	180	330

Figure 3. T7 Endonuclease I assay result for control

Lane M: 100bp plus DNA ladder (#M01010A).

Lane 1: unpurified PCR product from control, undigested

Lane 2: unpurified PCR product from control, digested

4. Search for potential SNP sites using online tools

Mammalian cells contain at least 2 copies of every chromosome, one copy from the mother, and the other from the father. So although most stretches of the chromosome will be identical to one another, it is possible to have some variations between them-in the forms of SNPs or other base pair differences in your target region sequence. Also, some tumor genes, such as P53, tend to have lots of mutations. So even in the negative control, it is possible to get T7 cutting from PCR product denaturation and reannealing. This is why when designing target primers, your design strategy should include avoiding SNPs in the first place.

We highly recommend you obtain sequence information of the target site of your cell line before primer design and any other experiment. NCBI provides a [tutorial](#) for searching for SNP information of a gene. We offer a brief description of the process here.

●BY GENE NAME

1. Search the [Gene](#) database with the gene name. If you know the gene symbol and species, enter them as follows: tpo[sym] AND human[orgn]



NCBI Resources How To Sign In


Gene

Advanced

Filters activated: Current only. [Clear all](#)

Gene

2. Click on the desired gene.

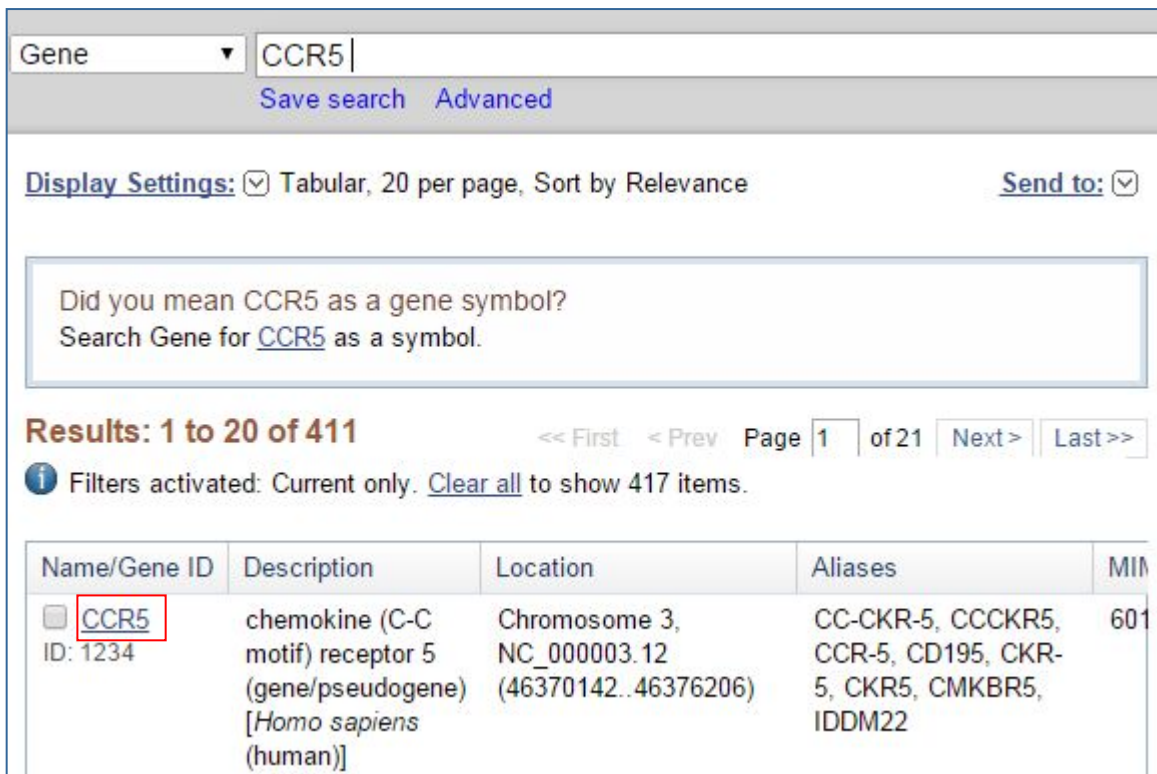

Gene

[Display Settings:](#) ☒ Tabular, 20 per page, Sort by Relevance [Send to:](#) ☐

Did you mean CCR5 as a gene symbol?
Search Gene for [CCR5](#) as a symbol.

Results: 1 to 20 of 411 << First < Prev Page 1 of 21 Next > Last >>

Filters activated: Current only. [Clear all](#) to show 417 items.

Name/Gene ID	Description	Location	Aliases	MIM
<input type="checkbox"/> CCR5 ID: 1234	chemokine (C-C motif) receptor 5 (gene/pseudogene) [Homo sapiens (human)]	Chromosome 3, NC_000003.12 (46370142..46376206)	CC-CKR-5, CCCKR5, CCR-5, CD195, CKR-5, CKR5, CMKBR5, IDDM22	601

3. In the list of Links on the right, click "GeneView in dbSNP". If the link is not present, no SNPs are currently linked to this gene.
4. For human genes, another option is to go to the variation section (Click on Variation in the table of contents in the upper right), and follow links to [Variation Viewer](#) for either the GRCh37/hg19 or GRCh39/h38 assemblies, to the [1000 Genomes Browser](#), [ClinVar](#) and more.

CCR5 chemokine (C-C motif) receptor 5 (gene/pseudogene) [*Homo sapiens* (human)]
Gene ID: 1234, updated on 17-May-2015

Summary

Official Symbol CCR5 provided by [HGNC](#)

Official Full Name chemokine (C-C motif) receptor 5 (gene/pseudogene) provided by [HGNC](#)

Primary source [HGNC:HGNC:1606](#)

See related [HPRD:03223](#); [MIM:601373](#)

Gene type protein coding

RefSeq status REVIEWED

Organism [Homo sapiens](#)

Lineage Eukaryota; Metazoa; Chordata; Craniata; Vertebrata; Euteleostomi; Mammalia; Eutheria; Euarchontoglires; Primates; Haplorhini; Catarrhini; Hominidae; Homo

Also known as CCR5; CCR-5; CD195; CKR-5; CCCR5; CMKBR5; IDDM22; CC-CKR-5

Summary This gene encodes a member of the beta chemokine receptor family, which is predicted to be a seven transmembrane protein similar to G protein coupled receptors. This protein is expressed

Table of contents

[Summary](#)

[Genomic context](#)

[Genomic regions, transcripts, and products](#)

[Bibliography](#)

[Phenotypes](#)

[Variation](#)

[HIV-1 interactions](#)

[Pathways from BioSystems](#)


[Interactions](#)

[General gene information](#)

[Markers, Clone Names, Homology, Gene Ontology](#)

[General protein information](#)

[NCBI Reference Sequences \(RefSeq\)](#)

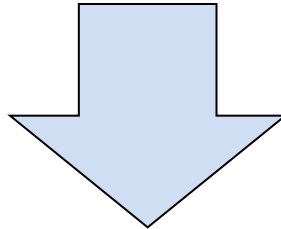

Variation

[See variants in ClinVar](#)


[See studies and variants in dbVar](#)

[See Variation Viewer \(GRCh37.p13\)](#)

☐ Genotypes

[See SNP Geneview Report](#)

[See 1000 Genomes Browser \(GRCh37.p13\)](#)

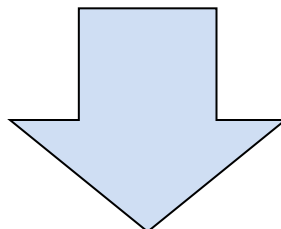

SNP linked to Gene (geneID:1234) Via Contig Annotation

The SNP GeneView page only reports human variation on GRCh38. A new [Variation Viewer](#) is available to view the gene CCR5 variations in [GRCh37p13](#) or [GRCh38](#), and will replace SNP GeneView later this year. Please visit the [Help Page](#) or [YouTube](#) for available features and send your comments and suggestions to NCBI [helpdesk](#).

Send

rs# on all gene models to Batch Query

Download

all rs# to file.

Gene Model (mRNA alignment) information from genome sequence

Total gene model (contig mRNA transcript):

2

mRNA	transcript	protein	mRNA orientation	Contig	Contig Label	List SNP
NM_000579.3	plus strand	NP_000570.1	forward	NT_022517.19	GRCh38.p2	<- currently shown
NM_001100168.1	plus strand	NP_001093638.1	forward	NT_022517.19	GRCh38.p2	View snp on GeneModel

☐ Clinical Source

☐ in gene region

☒ cSNP

☐ has frequency

☐ double hit

refresh

gene model	Contig Label	Contig	mRNA	protein	mRNA orientation	transcript	snp count
(contig mRNA transcript):	GRCh38.p2	NT_022517.19	NM_000579.3	NP_000570.1	forward	plus strand	259, coding

Region	Chr. position	mRNA pos	dbSNP rs# cluster id	Heterozygosity	Validation	MAF	Allele origin	3D	Linkout	Function	dbSNP allele	Protein residue	Codon pos	Amino acid pos	PubMed
	46372912	367	rs748244565	N.D.						missense	G	Glu [E]	1	4	
										contig reference	C	Gln [Q]	1	4	
	46372913	368	rs745912425	0.000						frame shift	-	Arg [R]	2	4	
										contig reference	AA	Gln [Q]	2	4	
	46372915	370	rs763192695	0.000						missense	A	Met [M]	1	5	
										contig reference	G	Val [V]	1	5	
	46372916	371	rs766432600	0.000						missense	G	Gly [G]	2	5	
										contig reference	T	Val [V]	2	5	
	46372934	389	rs751603911	0.000						missense	G	Gly [G]	2	11	
										contig reference	A	Asp [D]	2	11	
	46372936	391	rs755442066	0.000						missense	G	Val [V]	1	12	
										contig reference	A	Ile [I]	1	12	
	46372937	392	rs781613165	0.000						missense	C	Thr [T]	2	12	
										contig reference	T	Ile [I]	2	12	
	46372939	394	rs753095965	0.000						missense	G	Asp [D]	1	13	

Figure 4. Example of SNP Geneview Report on NCBI. Check column "mRNA pos" for your SNP site of interest. Check column "Chr. position" or "db SNP rs# cluster id" for the sequence of the site of interest.

VIII. Limited Use License and Warranty

Limited use license

Following terms and conditions apply to use of the IndelCheck™ CRISPR/TALEN insertion or deletion detection system. If the terms and conditions are not acceptable, the Product in its entirety must be returned to GeneCopoeia within 5 calendar days. A limited End-User license is granted to the purchaser of the Product. The Product shall be used by the purchaser for internal research purposes only. The Product is expressly not designed, intended, or warranted for use in humans or for therapeutic or diagnostic use. The Product must not be resold, repackaged or modified for resale, or used to manufacture commercial products or deliver information obtained in service without prior written consent from GeneCopoeia. Use of any part of the Product constitutes acceptance of the above terms.

Limited warranty

GeneCopoeia warrants that the Product meets the specifications described in the accompanying Product Datasheet. If it is proven to the satisfaction of GeneCopoeia that the Product fails to meet these specifications, GeneCopoeia will replace the Product. In the event a replacement cannot be provided, GeneCopoeia will provide the purchaser with a refund. This limited warranty shall not extend to anyone other than the original purchaser of the Product. Notice of nonconforming products must be made to GeneCopoeia within 30 days of receipt of the Product. GeneCopoeia's liability is expressly limited to replacement of Product or a refund limited to the actual purchase price. GeneCopoeia's liability does not extend to any damages arising from use or improper use of the Product, or losses associated with the use of additional materials or reagents. This limited warranty is the sole and exclusive warranty. GeneCopoeia does not provide any other warranties of any kind, expressed or implied, including the merchantability or fitness of the Product for a particular purpose.

GeneCopoeia is committed to providing our customers with high-quality products. If you should have any questions or concerns about any GeneCopoeia products, please contact us at 301-762-0888.

© 2015 GeneCopoeia, Inc.